

C-Series Syringe Pumps - 30 mm Stroke



PRODUCT OVERVIEW

- > OEM liquid handling component, delivers fluids from microliters to milliliters with exceptional accuracy and precision
- > Utilizes industry standard communication protocols and mounting configurations
- > C24000 is the preferred choice for low flow rate and low pulsation applications
- > Designed for continuous duty and reliable performance

FEATURES

- > Accommodates valves with 3 or 4 ports in a variety of configurations
- > High Resolution modes available, up to 192,000 microsteps/full-stroke
- > Uses industry standard OEM or DT Cavo® protocol
- > Communication interfaces: RS232, RS485 or CAN
- > Address up to 15 pumps individually
- > Fast valve switch time: < 250 ms between neighboring ports
- > Long-life rotary shear valves (up to 24 times the life of competing valves)
- > Proprietary high-efficiency and maintenance-free stepper motor drive system (no belts to stretch or break)

CONFIGURED OPTIONS

Fluid Contact Syringe Material	Borosilicate Glass and UHMWPE or PTFE (seal) - Pressures to 80 PSI
Fluid Contact Valve Material	PCTFE or PEEK (housing) and PTFE or UHMWPE (plug) - Pressures to 100 PSI
Syringe Sizes	50 µL, 100 µL, 250 µL, 500 µL, 1 mL, 2.5 mL, 5 mL, 12.5 mL

BASE MODEL

C3000
C24000




(Shown with accessory Syringe)

PRODUCT SPECIFICATIONS


	C3000 Syringe Pump	C24000 Syringe Pump
Stroke Speed	1.2 seconds to 100 minutes per stroke (800 minutes in microstep mode).	4.2 seconds to 800 minutes per stroke (3200 minutes in microstep mode).
Step Resolution	Standard mode 3,000 steps, Hi-Res mode 24,000 steps	Standard mode 24,000 steps, Hi-Res mode 192,000 steps
Drive Design	Direct drive stepper motor/leadscrew with 30 mm of travel. Manually moveable for "power off" syringe removal.	Direct drive stepper motor/leadscrew with 30 mm of travel.
Valve Options	120° 3-port, 90° 4-port, T-valve, 90° 3-way Distribution. Body materials PEEK or PCTFE, Rotary Plug PTFE or UHMWPE	
Interface Type	RS232, RS485 or CAN	
Communication Rates	RS232/RS485 9600 or 38400 baud CAN - 100K, 125K, 250K, 500K, 1M bps	
Format (RS232/RS485)	Data Bits - 8, Parity - None, Stop Bit - 1 Half Duplex, Data Terminal and OEM protocols	
Operating Temperature and Humidity	15°C to 40°C (59°F to 104°F) and 20% to 95% RH at 40°C (104°F), non-condensing	
Non-Operating Temperature and Humidity	-20°C to 65°C (-4°F to 149°F) and 30% to 85% RH, non-condensing	
Media Temperature	15°C to 40°C (59°F to 104°F)	

FUNCTIONAL PARAMETERS




Imprecision (full stroke) at 50 µL to 100 µL	< 0.1% CV
Imprecision (full stroke) at 250 µL to 12.5 mL	< 0.05% CV
Inaccuracy (full stroke)	< 1% syringe volume
Power Requirement	24 VDC ± 10%
Power Rating	≤ 1.5 Amp (peak)

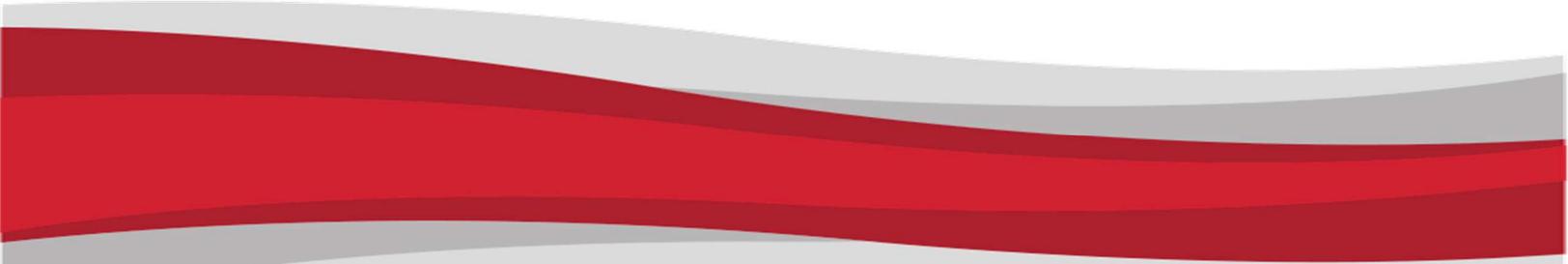

C-SERIES PUMP INTERFACE DETAILS

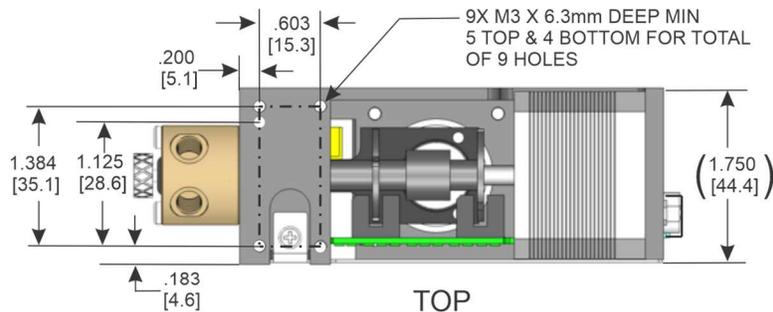

ALL MOUNTING HOLES ARE
M3 X 0.5 mm X 6.3 mm DEEP MIN

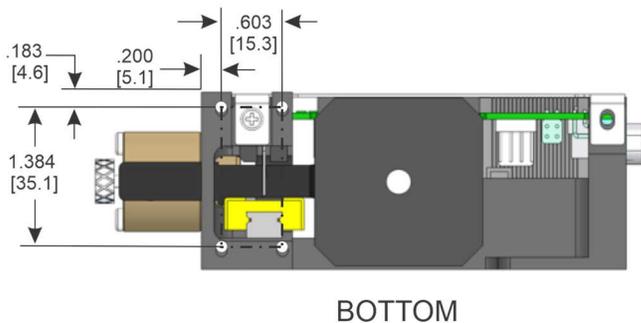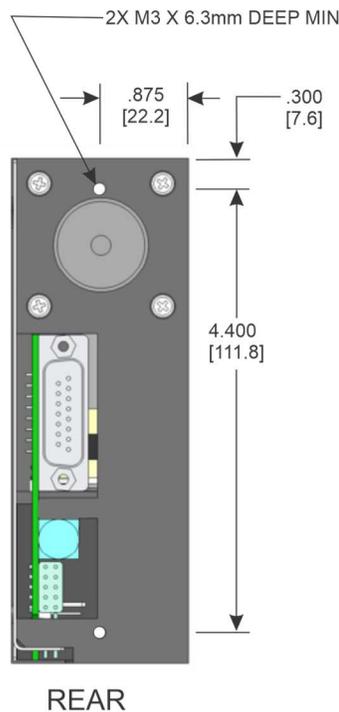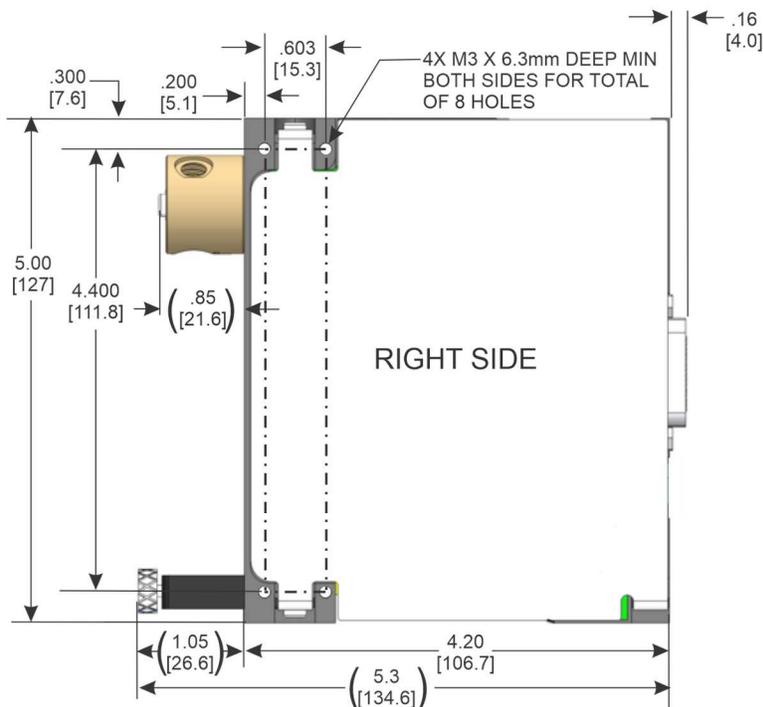

ALL DIMENSIONS UNLESS
OTHERWISE SPECIFIED
.XXX±.008 [0.2 mm]
XX±.016 [0.4 mm]
(MILLIMETERS SHOWN IN BRACKETS)
DIMENSIONS SHOWN IN PARENTHESES
FOR REFERENCE ONLY

PUMP AND VALVE SOLUTIONS FOR OEMS WORLDWIDE



tricontinent.com



Tricontinent Scientific, Inc.
12740 Earhart Ave
Auburn, CA 95602
USA
Tel: +1 800 937 4738
Fax: +1 530 273 2586
liquidhandling.tcs@gardnerdenver.com

Gardner Denver Thomas GmbH
Livry-Gargan-Str. 10
82256 Fürstenfeldbruck
Germany
Tel: +49 8141 2280 0
Fax: +49 8141 8892136
thomas.de@gardnerdenver.com

Gardner Denver Thomas Pneumatic Systems, Co., Ltd.
No. 1 New Dong An Road, Shuofang Town
Wuxi, Xinqu District, Jiangsu 214142
China
Tel: +86 510 6878 2258
Fax: +86 510 6878 2200
tricontinent.cn@gardnerdenver.com

The information presented in this material is based on technical data and test results of nominal units. It is believed to be accurate and reliable and is offered as an aid to help in the selection of Tricontinent products. It is the responsibility of the user to determine the suitability of the product for the intended use and the user assumes all risk and liability in connection there with. Tricontinent does not warrant, guarantee or assume any obligation or liability in connection with this information. Photos of products pictured in this catalog do not necessarily represent a specific model number. To obtain further information for custom options, contact your local Tricontinent office.

Printed in USA Form No. MKT90021 P © Gardner Denver Tricontinent. All rights reserved.